



# AVALO<sup>®</sup> SERIES AC MEDICATION CART

  
CAPSAHEALTHCARE

Avalo<sup>®</sup> ACS Medication Cart  
Avalo<sup>®</sup> ACM Medication Cart  
Avalo<sup>®</sup> ACXL Medication Cart



# Avalo® AC Medication Cart

## Smart, Secure & Simple

Capsa Healthcare brings its superior design expertise to the hospital market with the Avalo AC Medication Cart. This innovative product line features the latest design and serves as a platform for existing and emerging medication control technology.

Capsa has redefined medication cart security with a broad choice of lock systems including keyless entry with automatic relocking, mag stripe card access, keyless narcotic storage and cart user audit.

Best of all, AC Medication Carts offer huge capacity and intelligent solutions for medication management and organization.

### AC Medication Cart Highlights & Options

- A** Internal pass-thru waste receptacle
- B** Optional expansive slide-out writing surface
- C** Magnetic stripe, proximity or bar code security for cart access
- D** Optional keyless access with auto-relocking
- E** Integrated and flexible drawer divider system
- F** Configurable medication storage (2 bin sizes & 3 drawer sizes)
- G** Full selection of convenience accessories & options
- H** Stable, trim-line footprint for simple maneuverability
- I** Choice of sealed lead acid or LiFe battery system



## Cassette & Bin Options


**Designed with flexibility in mind,** Capsa Healthcare's cassettes and bins give you a choice of 2 or 3 tier durable cassettes with the option of 5.5" or 8" bins.

### 2-Tier & 3-Tier Cassettes



*The AC Medication Cart accommodates 2 & 3 tier cassettes to customize the storage configuration that best suits your requirements.*

### 8" & 5.5" Bins

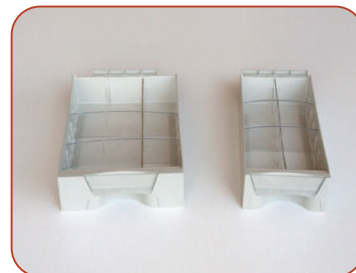

*AC cassettes give you an option of 8" and 5.5" bins that offer a flexible divider system. Large bin lens permits protected and clear patient labeling.*





# Avalo® AC Medication Cart

## Security Options

Capsa Healthcare has fitted AC Medication Carts with its proven and innovative automatic relocking system and the latest integrated medication cart security.

To further optimize security, add a secondary access requirement such as ID badge card, bar code or proximity reader.

### Lock System Choices

#### AutoLock™

Automatically secures a motionless unlocked medication cart with a user-defined relock time delay from 1-99 minutes. AutoLock permits up to 5,000 individual user codes to further enhance cart access tracking and security. *See available security upgrades.*



#### Keyless Access with Manual Relocking

Simply enter a 4-digit access code, push down the lock lever and the cart is ready for access. To secure the medication cart, simply push in the lock lever. Clear visual indication of lock status is apparent from a distance.



#### Key Access with Manual Relocking

Featuring the Best® Core Removable Lock System, this basic yet proven lock system provides your medication cart advanced security and the flexibility to change keys simply and easily.



### AutoLock™ Security Upgrades

#### Keyless Narcotics Storage



*The optional keyless narcotics storage feature greatly augments security to further optimize control and accountability.*

#### Secondary Access



*Select from a bar code, mag stripe or proximity card reader to serve as the primary or secondary access requirement*

#### Audit Tracking



*Authorized users can track cart activity, including user access and relocking for over 2,000 events.*

## Accessories & Options

### Customize your configuration


with AC Medication Carts. Streamlined clinical workflow is assured with a full complement of convenience accessories and organizing supplies to facilitate patient care.



IV Pole



Drawer Dividers

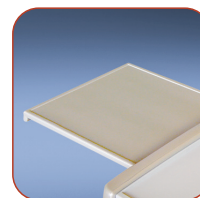

Writing Surface




Sharps / Gloves



Clear Top Mat



Cart Light

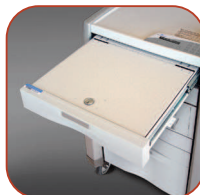

Narcotics Storage



Waste Container



Bulk Storage



Mobile Computing

# Avalo® AC Medication Cart

## Specifications & Models

### AC Medication Cart Features & Specifications

The Avalo AC Medication Cart line combines the latest medication cart design with plenty of capacity and simple maneuverability. Choose the drawer configuration and cassette size that best meets your requirements for highly-secure and organized medication control.

This solution is perfect for point-of-care medication administration and as a delivery workstation to augment your dispensing cabinet.

### AC Cart Models

Choose from 3 cart models. ACS and ACM available in standard or compact height, ACXL available in standard height only.



### Features

- Power status display on LCD
- Large retractable writing surface (optional)
- Internal pass-thru waste receptacle (ACS/ACM)
- Modular bin / storage drawers
- Expansive storage drawers
- Slot-Lock™ security system secures all partially closed drawers and bins left open up to 1 inch
- Keyless access and auto-relocking with AutoLock option
- 5" double ball bearing casters (12.7cm)

### Benefits

- Improved workflow and accuracy
- Safe and ergonomic design
- Advanced medication security and control
- Lightweight design for easy maneuverability
- Customized medication / supplies storage

### Standard Cart Specifications

<b>Base Dimensions:</b>	ACS: 24"d x 26"w / 61cm x 66cm ACM: 24"d x 38"w / 61cm x 96.5cm ACXL: 24"d x 48.5"w / 61cm x 123.2cm
<b>Work Surface Area:</b>	ACS: 16"d x 18"w / 40.5cm x 46cm ACM: 16"d x 28"w / 40.5cm x 71cm ACXL: 19.75"d x 41"w / 104cm x 50.1cm
<b>Work Surface Height:</b>	Compact: 39.5" / 100cm Standard: 43" / 109cm
<b>Construction:</b>	Formed steel, aluminum, heavy-duty polymer work surface / end panels
<b>Finish:</b>	High-density polymer
<b>Wheels:</b>	Swivel casters with 1 locking; tracking caster optional

Specifications, availability & components are subject to change without notice.

