

Tryten P5 Monitor Cart



DESIGNED FOR:


TELEHEALTH COMMUNICATION

REMOTE INTERPRETING

VIRTUAL VISITING

VIRTUAL ROUNDING

VIRTUAL PATIENT OBSERVATION

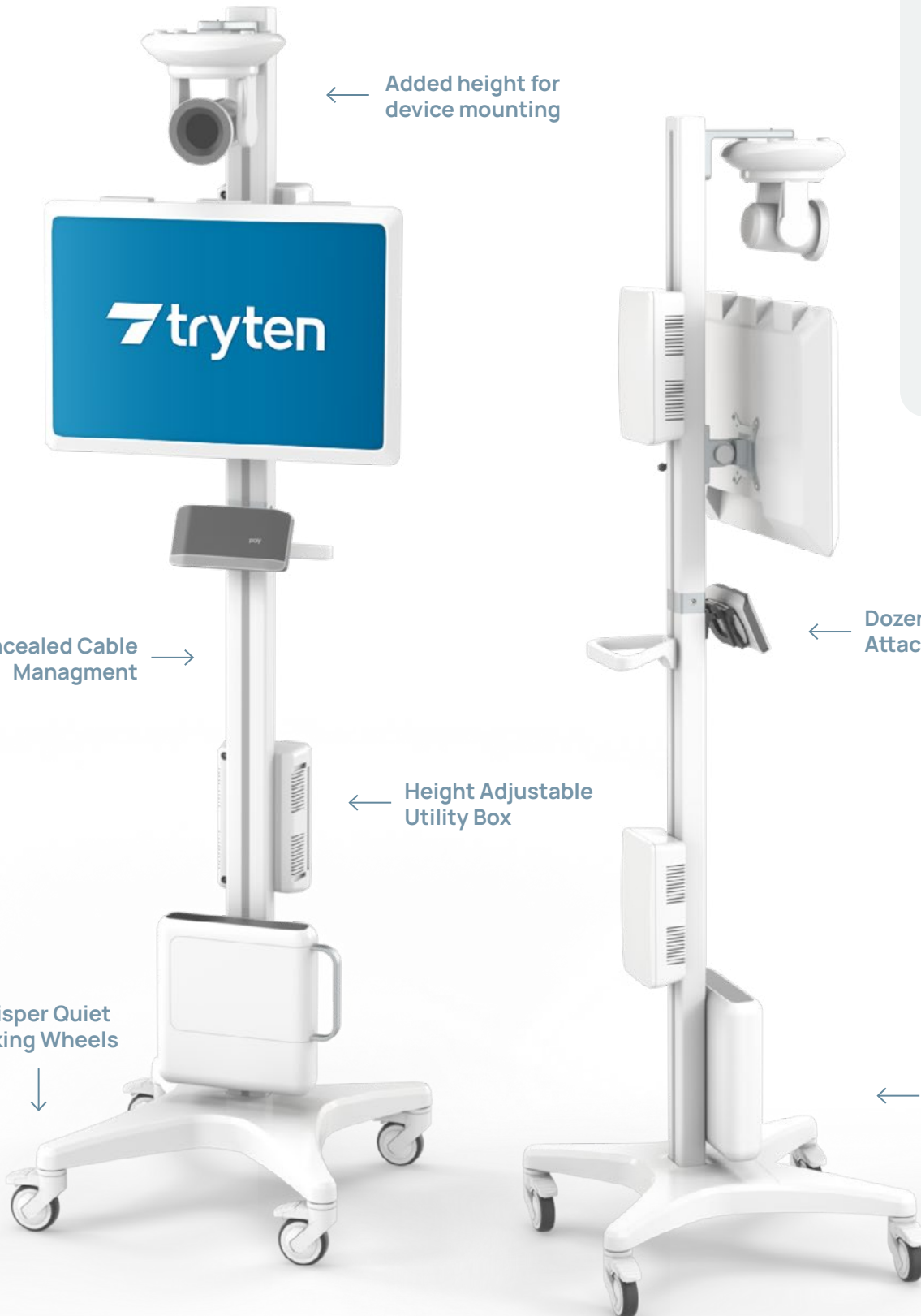

← Added height for device mounting

→ Concealed Cable Management

↓ Whisper Quiet Locking Wheels

← Height Adjustable Utility Box

← Dozens of Accessory Attachment Options

← Tip Resistant Weighted Base

Tryten P5 Monitor Cart



Wide Straight P-Series Pole

The new, larger, pole allows for the mounting of larger monitors, accessories and other peripherals.

Ergonomic Handle & Cord Wrap.

Ergonomic handle with slender design uniquely placed for easy push-and-position.

Best-in-Class, Durable Materials

Anodized Aluminum, Powder-Coated Aluminum, Cycloy Plastic, fully wipeable and approved in operating rooms.

Customization

Customize your cart with Branding Options and 40+ Accessories.

Concealed Cable Management

A wider column allows for cable management on both sides of the internal channel, for discreet and uninterrupted charging. Large utility box can house power cable receptacles or other power components.

Base & Casters

The 4-Wheel, larger base maximizes stability and balances the need for stability with a compact footprint. Complete with premium, quiet, non-marking, locking casters, this cart is ready to roll or stay in a locked position.

Quality & Services

5 Year Limited-Warranty, ISO 9001:2015 Certified, Tip Resistant as per IEC 60601, and Remote Expert Support.

DIMENSIONS AND SPECS

- Total cart height: 71.7"
- Column extension: 28"
- Large Utility Box, Inside: 8.6" x 3.9" x 2.0"
- Large Utility Box, outside, 9.2" x 5.5" x 2.8"
- Base footprint: 18.7" to 20.6"
- Maximum mounting weight 30lbs

WHAT'S INCLUDED

- P-Series Pole
- Straight Adapter and Extension
- Base with Locking Casters
- Large Utility Box
- Ergonomic Handle
- Hardware Pack
- Assembly Instructions