

Tryten P6 Monitor Cart




DESIGNED FOR:

TELEHEALTH COMMUNICATION

TELESTROKE

VIRTUAL ROUNDING

CONFERENCE ROOM ON WHEELS

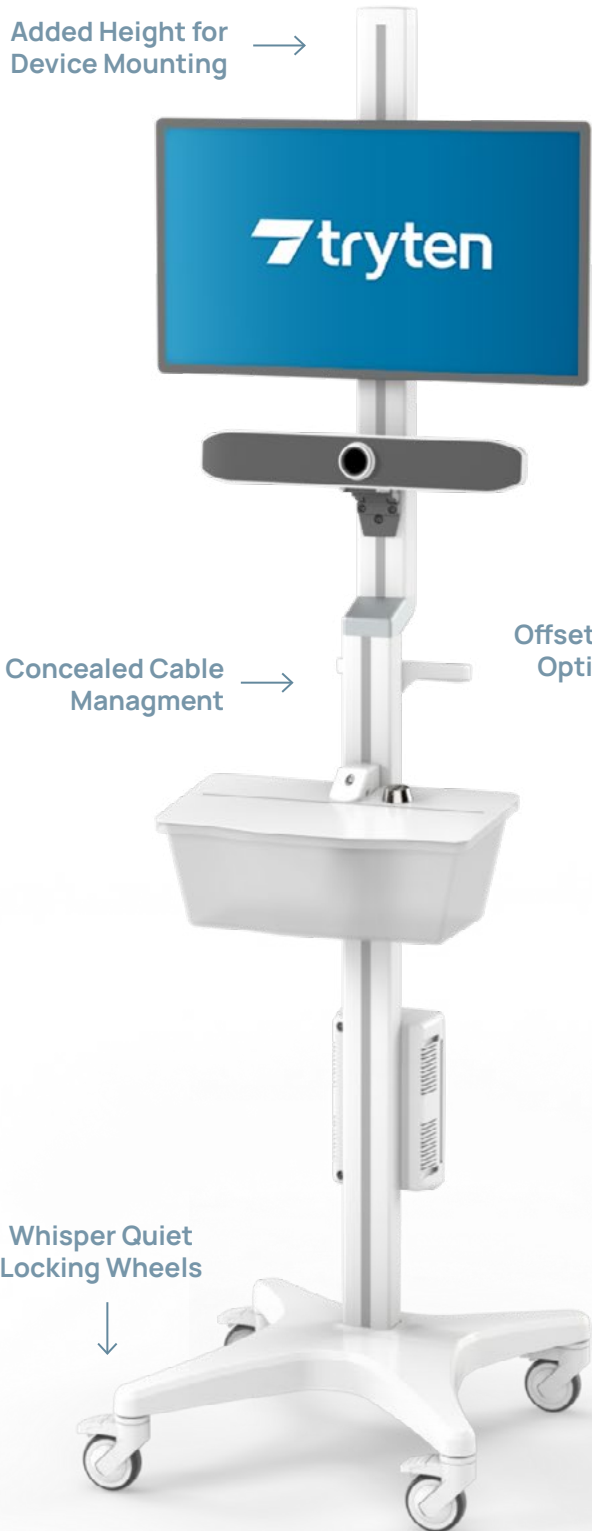

Offset Adaptor for Optimal Viewing



Tryten P6 Monitor Cart



Wide Offset P-Series Pole

The new, larger, pole allows for the mounting of larger monitors, accessories and other peripherals.

Ergonomic Handle & Cord Wrap

Ergonomic handle with slender design uniquely placed for easy push-and-position.

Best-in-Class, Durable Materials

Anodized Aluminum, Powder-Coated Aluminum, Cycloy Plastic, fully wipeable and approved in operating rooms.

Customization

Customize your cart with Branding Options and 40+ Accessories.

Concealed Cable Management

A wider column allows for cable management on both sides of the internal channel, for discreet and uninterrupted charging . Large utility box can house power cable receptacles or other power components.

Base & Casters

The 4-Wheel, larger base maximizes stability and balances the need for stability with a compact footprint. Complete with premium, quiet, non-marking, locking casters, this cart is ready to roll or stay in a locked position.

Quality & Services

5 Year Limited-Warranty, ISO 9001:2015 Certified, Tip Resistant as per IEC 60601, and Remote Expert Support.

DIMENSIONS AND SPECS

- Total cart height: 72"
- Column extension: 28"
- Large Utility Box, Inside: 8.6" x 3.9" x 2.0"
- Large Utility Box, outside, 9.2" x 5.5" x 2.8"
- Base footprint: 18.7" to 20.6"
- Maximum monitor mounting weight: 25lbs

WHAT'S INCLUDED

- P-Series Pole
- Offset Adapter and Extension
- Base with Locking Casters
- Large Utility Box
- Ergonomic Handle
- Hardware Pack
- Assembly Instructions